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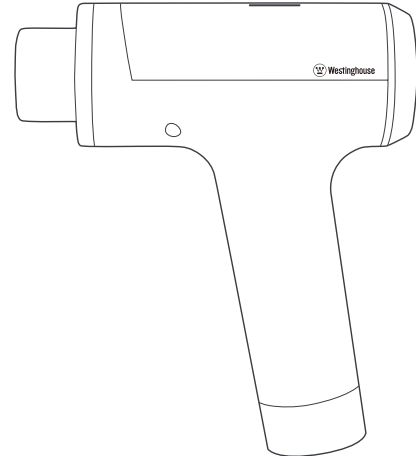
# Warnings

Before undertaking massage please first obtain approval from your doctor if any of the following apply: Pregnancy, diabetes with complications such as neuropathy or retinal damage, peripheral neuropathy, wearing of pacemakers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylolysis, recent joint replacements or IUDs, metal pins or plates or any other concerns about your physical health. Not for use by persons 13-year-old and under.

**Do not use this device on your head or near your genitals.** Do not use the device above your C4 (approximately mid-throat level).

**Do not place your finger or any objects near the shaft or socket while the device is in use.** Do not get hair caught in the massage head socket or shaft or any other moving part on the device.

**Do not block the ventilation ports of the motor or battery.** Do not allow dust or debris to build up in vent ports. Do not insert objects in the motor or battery ventilation ports.

**Do not get the device wet other than by cleaning it with a lightly damp towel.**

**This device contains an extremely powerful magnet!** Users of pacemakers must consult their doctor before using this device.

**Do not use any batteries other than the ones provided by Flow Sports Technology.** Only use as directed by these instructions.

# Safety


- Those aged 14 years and above or persons with reduced physical, sensory or reasoning capabilities should only use this device once they have received supervision or instruction by a responsible person who has demonstrated use of the device in a safe way and disclosed the hazards involved.
- Do not use device as a toy. Close attention is necessary when used near children.
- Keep the device out of reach of children.
- Never submerge the device in water or other liquids or pour water over the device.
- Keep away from heat sources and direct sunlight.
- Keep ventilation ports free from dirt, dust, and debris.
- Do not use the device for any other purpose than personal sports massage.
- Do not handle any part of the A/C wall charger or device with wet hands.
- Do not unplug the A/C wall charger by pulling on the cable. Grasp the plug to unplug, not the cable.
- Do not use any lubricants, cleaning agents, polishes or air fresheners on any part of the device.
- Contact the Flow support team when repair is required. Do not disassemble the device as incorrect reassembly may result in an electric shock or fire.
- If the device is malfunctioning, or has been dropped, damaged, left outdoors, received a sharp blow or dropped into water, do not use and contact the Flow support team.

# Battery Safety

- Do not connect battery poles with conductors such as keys, coins, jewellery, or other metal materials. Doing so may short-circuit the battery causing injuries.
- Keep the battery away from excessive heat and direct sunlight. Do not place battery on or in heating devices, such as microwave ovens, stoves, or radiators. Batteries may explode if overheated.
- Do not put battery in a fire as it may explode. Damaged batteries may also explode. Dispose of used batteries in accordance with local regulations.
- Improper battery use may lead to a fire, explosion, or other hazards.
- Do not allow children or pets to bite or suck the battery. Doing so may result in damage or explosion.
- Do not smash or pierce the battery, or expose it to high external pressure. Doing so may lead to a short circuit or overheating.
- Replacing the battery with an incorrect model or non-genuine copy may cause an explosion.

# In the Box


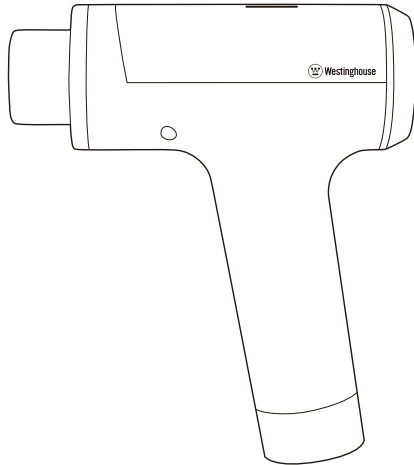
Also included:

- Soft carry bag
- Operating guide


A/C Wall Charger

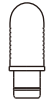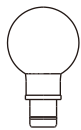
Heads

# Overview
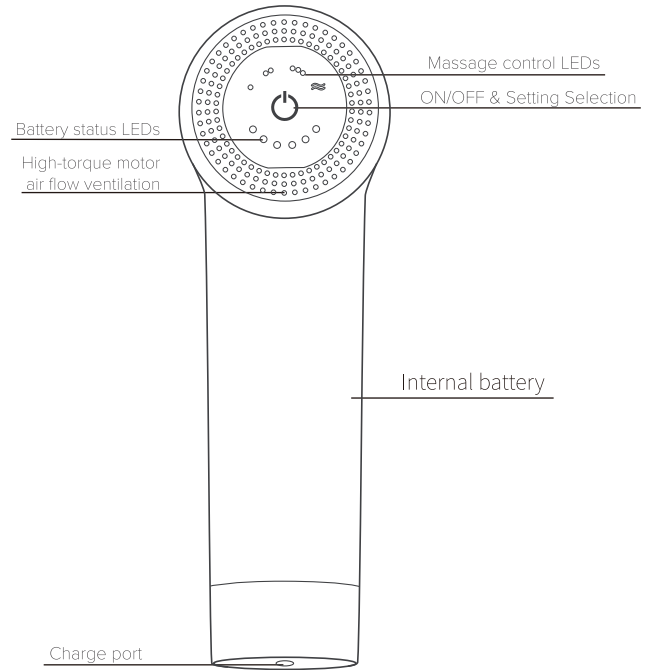
# Heads

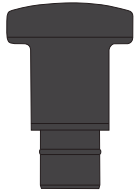
### Hammer

The jack of all trades head, this head is designed for overall body treatment and larger muscles groups.


### Bullet

The small diameter of the bullet head makes it perfect for focused muscle treatment, to pinpoint trigger points and work them away.
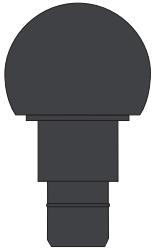
### Medium


Made with closed-cell PU foam, this forgiving, non-porous & hygienic head is perfect for all over body use and more sensitive zones.

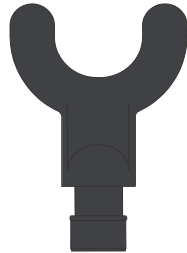
### Fork

Perfect for the spine, lower neck and Achilles' tendon along with major muscle groups like the calves and quadriceps. Care should be taken when using the fork head in tender areas.

# Quick Start

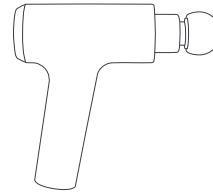
Remove the device carefully from the box. Device will already have the Medium head fitted. To remove it simply grip the base of the head and tug until you hear a 'pop'.

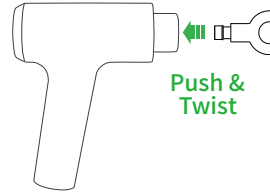**Push & Twist**


With the device off, fully insert the Head of your choice by pushing it into the socket. Do not operate the device without a head.

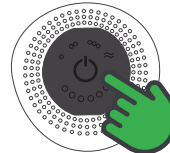**Press**


To turn the device on press the large power button in the centre of the control panel. Pressing again will activate the various modes (see opposite). Massage desired body part at the required speed and apply pressure as needed. Cycle through the modes to switch off when finished.

# Modes

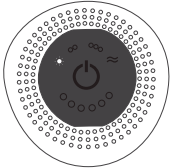
Your device is designed to be as easy to operate as possible.


Press the large Power button once, and the device will start to operate in **Relax mode**, our first level.

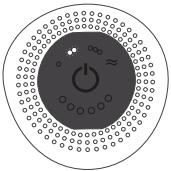Press again for level two, **Refresh mode**.

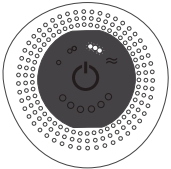Press once more for maximum power at level 3, **Recharge mode**.

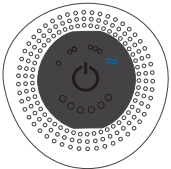Press a fourth time to start **Conversion mode**.


To switch off press one last time and your device will power down. There is no need to do anything else, as it will sleep after 5 minutes of non-use. To wake, simply press the Power button.

# Charging

## Low or No Charge

When your device has very low or zero charge left in its battery the control panel will display a **Red** LED and may operate sluggishly.


It is not recommended to run the battery down to this point.
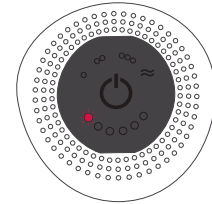
## Charging your Device

Plug the A/C wall charger securely into a wall socket.


Press power button to charge device. The charger unit will display a **Red** LED indicating that the charger is connected and providing charge - the charge level LEDs on the control panel will cycle through - **Red, Amber & Green**.
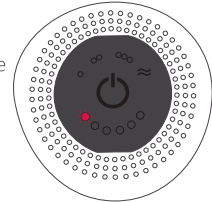
## Charging Complete

When fully charged the device will show all 6 power indicator LEDs filled.

The device will not turn on while charging.
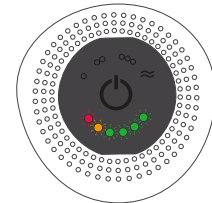
# Cleaning

- Do not clean device while charging.
- Do not use water to clean the device.
- Do not use abrasive cleaning chemicals or solvents.
- Do use a slightly damp towel to clean the device.
- Do use the supplied cloth to dry afterward.

# Maintenance

- High-density foam heads should be replaced every three months or if they show obvious signs of wear and tear or become dirty.
- Fully charge battery before first use.
- Recharge device monthly to extend battery life.
- Device can be recharged at any time and at any battery level.
- Only recharge using the supplied A/C Wall Charger.
- Do not keep/leave/store the device in direct sunlight or at temperatures above 40°C.
- The device and charger work best at normal room temperature. Do not charge in direct sunlight or cold temperatures. The temperature range for charging is 0°C to 40°C.
- Never leave your device operating or charging unattended.
- It is recommended not to fully deplete the battery.

# Travel/Storage

- It is TSA approved as a carry-on item. It must only be stored in carry-on baggage. If a carry-on bag is checked at the gate, all spare lithium batteries must be removed from the bag and kept with the passenger in the aircraft cabin.
- Fully charge the battery prior to long-term storage.

# Warranty

## TERMS AND CONDITIONS OF THE 12 MONTH LIMITED WARRANTY

### WHAT IS COVERED

- The repair or replacement if it is found to be defective due to faulty materials, workmanship or function within 12 months of purchase or delivery.
- Where this device is sold outside of the EU, this warranty will only be valid if the device is used in the country in which it was sold.
- Where this device is sold outside of the EU, this warranty will only be valid if the device is used in the country in which it was sold.

### WHAT IS NOT COVERED

Not guarantee the repair or replacement of a product where a defect is a result of:

- Accidental damage or faults caused by negligent use or care, misuse, abuse, neglect, carelessness using the device in a way which is not in accordance with the Operating Manual.
- Damage caused by not carrying out the recommended device maintenance.
- Damage caused by the failure to adhere to device cleaning directions.
- Damage caused by the failure to adhere to travel and storage directions.
- Use of the device for anything other than as a handheld sports massager.
- Use of parts not assembled or installed in accordance with the instructions.
- Use of parts and accessories which are not genuine components.
- Normal wear and tear.

If you are in any doubt as to what is covered by your warranty, please contact the customer support team.

### PLEASE NOTE THIS WARRANTY REQUIRES THE DATED PROOF OF PURCHASE

If you are a reseller partner, please contact the appropriate wholesaler, dealer or retailer, from whom you directly purchased the product for warranty and return authorisation procedures.